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STERILE TRANSFUSION SET FOR TRANSFUSION OF BLOOD HEMOPLAST®, PK 21-02

Purpose

Sterile transfusion set for transfusion of blood HEMOPLAST®, PK 21-02 (transfusion set), hereinafter - set, is intended for transfusion to recipient of blood and blood components.

General information

Set is used in clinics, hospitals, first-aid service, field conditions.

Set is intended for transfusion from bottle.

Set is sterile, non-pyrogenic, non-toxic and intended only for single use.

Set works only under action of gravity.

20 drops of distilled water are equivalent to (1,0 ± 0,1) ml.

Connecting cone of set head - "Luer", which may have element of closing connection - "Luer-Lock".

Set may be supplied with injection needle.

Set may be supplied without injection unit.

Set may be completed with plaster for fixation of intravenous catheter.

ATTENTION!

DO NOT USE if integrity of consumer packaging is damaged and after expiry date.

DO NOT USE if one part of set is separated or protective cap is disconnected.

DO NOT USE under pressure more than 20 kPa after puncture of injection unit.

Contraindications

Do not use to people with natural latex allergy.

Side effects

It may cause allergic reaction in people who have increased sensitivity to natural latex.

It may increase risk of endocrine system malfunction in children and pregnant women with prolonged use.

Rules for use

Check serviceability of transfusion media for use.

Stir transfusion media in bottle. Cleanse plug of bottle by alcohol, iodine or any other antiseptic.

Take out set and put it on sterile table.

Close clamp 3. Take off cap 1 from needle 4 to connect to bottle. Introduce

needle 4 into plug of bottle up to stop. Turn over bottle and fix it on rack.

Fill in drip chamber 6 till immersion of filter into transfusion media, periodically pressing on lower part of drip chamber body.

Release excess pressure from bottle. For this take off cap 8 from head 5.

Open clamp 3.

Fill in transfusion media into tube of set, at this time head 5 should be at level of drip chamber 6.

Close clamp 3. Put cap 8 on head 5.

Take off cap 1 from needle 2 of air conduit. Introduce needle 2 into plug of bottle up to stop.

Take off cap 8 from head 5.

Open clamp 3 and slowly fill in set till full displacement of air from tube and injection unit 7 and emergence of drop from head 5.

Close clamp 3. Put cap 8 on head 5.

Check lack of air bubbles in tube of set.

When using set with injection needle, fulfill all above indicated actions with injection needle instead of cap 8.

When connecting set to catheter which is in vein of patient, cleanse plug of catheter by antiseptic. Take off plug from head of catheter and cleanse head by antiseptic. Drain solution of heparin from catheter until emergence of blood drop on head exit. Take off cap from the head "Luer" or "Luer-Lock" of set and connect it to head of catheter.

Use set for its intended purpose according to instructions of doctor.

Adjust required transfusion rate by means of clamp 3.

Introduce auxiliary medicines only through injection unit 7, using needle with diameter not more than 0,8 mm with closed clamp 3.

When changing bottle, close clamp 3 without withdrawal of injection needle or catheter from vein and not permitting discharge of drip chamber 6, take out needles 4 and 2 from empty bottle and introduce them quickly into plug of filled bottle - first needle 4, then needle 2.

After finishing of intravenous injections close clamp 3 and take out injection needle from vein of patient. When conducting intravenous injections through catheter, disconnect set from catheter, fill in catheter by solution of heparin and screw plug on head of catheter.

Guarantee of manufacturer

Expiry date - 5 years from date of manufacture.

Storage conditions

Store at temperature from 5 °C to 40 °C and relative humidity not more than 80 % at temperature of 25 °C.

Utilization

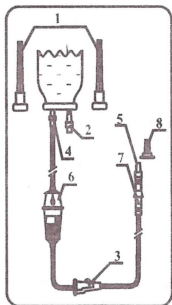
Utilize sets according to national requirements.

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